

# JoinPlant



 **MEDENCY**  
TECHNOLOGY | INNOVATION | PASSION


**JOINPLANT** is the fast and simple concept to fix implants in the event of immediate loading. Implant stability plays a critical role for successful osseointegration. This high-tech intraoral welding machine avoid implant mobility throughout the entire process of rehabilitation until osseointegration. Thanks to the exclusive CTP (Continuous Twin Pulse), **JOINPLANT** connect a titanium wire to the abutment. In few milliseconds, a high intensity electric charge is delivered to the contact point without any sparks or heating risk.

One touch pre-set selection and dedicated intraoral forceps make the procedure very quick and safe. **JOINPLANT** reduces the number of appointments with less treatment costs than any other technique. It works with all dental implants without any limitation

### GET TRAINED FROM EXPERIENCED SURGEON

Throughout academies and universities programs, Medency offers leading-edge training courses for doctors worldwide to help them successfully enhance their practice with integrative, proven and cost-effective techniques.

Joinplant Welding System	
Code	MEDMA001.1
Welding / Peak / Power	1 kW
Frequency	50 kHz
Energy	20 J - 200 J
Pulse range	Continuous Twin Pulse
Medical / Laser class	IIB
Dimension / Weight	230 X 200 X 170 mm / 6.0 kg

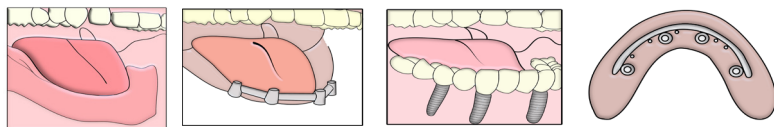

Immediate reinforcement of the temporary by a titanium framework. Safely and predictably for immediate loading treatment. Significant time and cost reduction.



For medical professional only

©2017 MEDENCY. All rights reserved. DPDENE04.0

